

Quality Built, American Made  
*Medical Equipment*

MADE  
IN THE **USA** 

**2026 Product Catalog**

[www.macmedical.com](http://www.macmedical.com)



## About Us

MAC Medical, Inc. manufactures top quality medical equipment including Warming Cabinets, Sinks, Tables, Casework, Case Carts, IV Stands, Mayo Stands, Linen Hampers, and many other stainless steel products.

To complement our stainless steel products, we have expanded our product line to include Plastic Laminate Casework and carts for a variety of spaces including medical, dental, laboratory, and more.

Over the last two and a half decades, MAC Medical has grown into a strong team dedicated to designing and manufacturing high-quality medical equipment. Located in Southern Illinois, our facilities provide an ideal hub for centralized distribution across the United States. Not limited to the United States, MAC Medical products can be found in over 20 countries worldwide.

Founded as a solutions-based company, MAC Medical continues to make that its core competency. With no project too large or too small, MAC Medical is committed to providing solutions for customer needs through superior service and American-made products.



### Corporate Office

609 S. Breese St.  
Millstadt, Illinois 62260

### Manufacturing Plant

820 South Mulberry Street  
Millstadt, Illinois 62260

Phone: 618-476-3550  
Toll Free: 877-828-9975  
Fax: 618-476-3337  
sales@macmedical.com  
www.macmedical.com



@macmedicalinc



A philosophy, a mindset, a mission statement...whatever you want to call it, our employees at MAC Medical are "Inspired To Be Better". It's this mentality that helps us to provide the quality built, American made products that our customers have come to expect.

Call Toll Free: 877-828-9975

## Table of Contents

<b>Warming Cabinets</b>	<b>2</b>	<b>Stainless Steel Tables</b>	<b>36</b>
Warming Cabinets Overview	3	Prep & Pack Workstations	37
Single Chamber Units	4	Utility Tables	38
Dual Chamber Units	5	Nested Tables	38
Triple Chamber Units	6	Work Tables	39
Warming Cabinet Accessories	6	Space Saver Tables	39
Warmer Options & How to Order	7	Adjustable Over Operating Tables	40
		Central Supply Tables	40
<b>Transport Stretchers</b>	<b>8</b>	Instrument Tables	41
General Transport Stretcher	9		
5th Wheel Transport Stretcher	10	<b>Stainless Steel Stands</b>	<b>42</b>
Stretcher Comparison	11	Mayo Stands <i>Redesigned</i>	43
		IV Stands	44
<b>Surgical Scrub Sinks</b>	<b>12</b>	Solution Stands	45
Surgical Scrub Sinks (ES-Series)	13		
Surgical Scrub Sinks (SS-Series)	14	<b>Other Stainless Steel Products</b>	<b>46</b>
Surgical Scrub Sinks (SS-Series) Options	15	Pass Through Windows	47
		Kick Buckets	47
<b>Processing Sinks</b>	<b>16</b>	Gowning Benches	48
Processing Sinks Overview	17	Step Stools	48
Processing Sinks Accessories	18	Linen Hampers	49
		Stationary Wall Shelves	50
<b>Stainless Steel Casework &amp; Cabinets</b>	<b>20</b>	Fold Down Wall Shelves	50
Wall Cabinets	22	Peg Boards	51
Base Cabinets	22		
High Cabinets	22	<b>Mobile Medical Carts</b>	<b>52</b>
Operating Room Cabinets	23	Short Mobile Carts	53
Casework Accessories	24	Tall Mobile Carts	54
Nurses Documentation Stations	26	Mobile Cart Accessories	55
Mobile Cabinets	27		
		<b>Green Initiative</b>	<b>56</b>
<b>Stainless Steel Carts</b>	<b>28</b>	<b>Warranty Information</b>	<b>57</b>
Closed Case Carts	29		
Sealed Case Carts	30	<b>Contact Us</b>	<b>58</b>
Vertical Handle Case Carts	31		
Open Case Carts	32		
Utility Carts	34		
Sterile Wrap Carts	35		




# Warming Cabinets

MAC Medical offers Warming Cabinets to help compliment and enhance your facility. As with all of our products, warming cabinets undergo stringent testing procedures during the manufacturing process to assure a long life, the highest quality, and accurate temperature settings. We offer standard sizes along with custom sizes to meet specific requirements.

Our Warming Cabinets provide independent, digitally controlled heating chambers that offer actual temperature and set point displays. Data is recorded every 30 minutes and stored for up to two years. Choose from single, dual, or triple chamber units.



**DWC243074E-G-4B**  
(shown with optional glass doors)

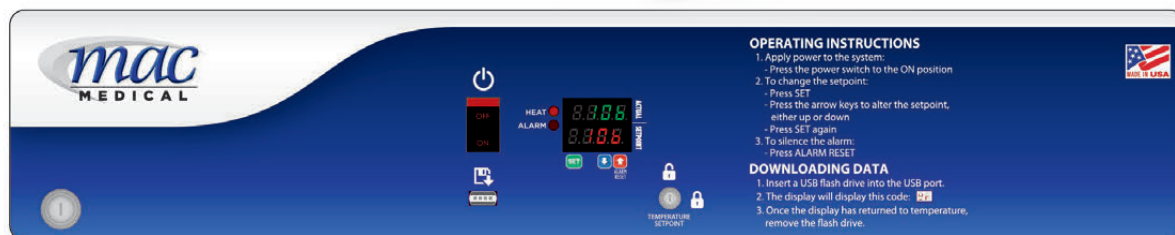

Header Control (Single Chamber)

## STANDARD FEATURES

- "Easy Log" integrated data logger with removable USB
- Keyed temperature lock out
- Keyed door lock for single & upper chamber units
- Adjustable feet
- Stainless steel doors
- Automatically records data every 30 minutes
- Internal storage up to 2 years
- Displays actual temperature and set points at all times
- Independent digitally controlled heating chambers
- Rapid warm time (2-6 hours depending on load)
- Push button digital operation
- Large digital display(s)
- Fully insulated to provide uniform heating
- Audible & visual over-temp alarm

## OPTIONAL FEATURES

- Full view glass door
- Left hand door hinge
- 220/240V power supply
- Celsius
- Two (2) recessed unit options
- Seismic braces
- Pass through chamber
- Electronic keypad lock
- Sloping top
- Keyed door lock for intermediate & lower chamber units
- Custom sizes available

## To Download Data

1. Insert supplied USB flash drive into USB slot.
2. Data download will automatically begin.
3. During download display will read:



4. Once the display shows the temperature again, it is safe to remove the flash drive.
5. Data is saved on the flash drive as a .CSV file for easy import into most spreadsheet programs (e.g. Microsoft Excel).
6. Data is secure (no reporting or connection to any network).

## Sample Report

SWC243074 Single, S/N: 062112-D01, 101118

DATE/TIME	TEMPERATURE	°C/°F
10/10/2018 23:12	106	F
10/10/2018 23:43	106	F
10/11/2018 0:12	106	F
10/11/2018 0:43	106	F
10/11/2018 1:12	106	F

# Single Chamber Units

## Warming Cabinets

### SWC151822

#### Overall Dimensions

17.25" D x 18" W x 22" H  
(438.15mm x 457.2mm x 558.8mm)  
Optional Adjustable Shelf

#### Chamber I.D.

14" D x 14" W x 11.63" H  
(355.6mm x 355.6mm x 259.4mm)  
1.30 Cubic Feet Storage  
Temperature Range 90°-160° F (Tolerance +/- 1° F)

#### Electrical

120V, 2.08 Amp, 50/60 Hz, .25 kWh (Avg.),  
853 BTU/hr (Avg.)



SWC151822-G-NB  
No base (NB) standard for this model

Additional items available for this model (sold separately):

- Mobile Stand (MS1518)
- Mobile Cabinet Stand (WBC1518)
- Roll Out Basket (WB1518)

### SWC182424

#### Overall Dimensions

20.5" D x 24" W x 24.5" H  
(520.7mm x 609.6mm x 622.3mm)  
Optional Adjustable Shelf

#### Chamber I.D.

17" D x 20" W x 15.25" H  
(431.8mm x 508mm x 387.3mm)  
3.0 Cubic Feet Storage  
Temperature Range 90°-160° F (Tolerance +/- 1° F)

#### Electrical

120V, 2.9 Amp, 50/60 Hz, .23 kWh (Avg.),  
785 BTU/hr (Avg.)



SWC182424-G-NB  
No base (NB) standard for this model

Additional items available for this model (sold separately):

- Mobile Stand (MS1824)
- Mobile Cabinet Stand (WBC1824)
- Roll Out Basket (WB1824)

### SWC243024 *Best Seller!*

#### Overall Dimensions

26.5" D x 30" W x 24.5" H  
(673.1mm x 762mm x 622.3mm)  
1 Adjustable Shelf

#### Chamber I.D.

23" D x 26" W x 15.25" H  
(584.2mm x 660.4mm x 387.3mm)  
5.27 Cubic Feet Storage  
Temperature Range 90°-160° F (Tolerance +/- 1° F)

#### Electrical

120V, 6.3 Amp, 50/60 Hz, .45 kWh (Avg.),  
1535 BTU/hr (Avg.)



SWC243024-G-NB  
No base (NB) standard for this model

Additional items available for this model (sold separately):

- Mobile Stand (MS2430)
- Mobile Cabinet Stand (WBC2430)
- Roll Out Basket (WB2430)

### SWC243036

#### Overall Dimensions

26.5" D x 30" W x 36.5" H  
(673.1mm x 762mm x 927mm)  
2 Adjustable Shelves

#### Chamber I.D.

23" D x 26" W x 25" H  
(584.2mm x 660.4mm x 635mm)  
8.65 Cubic Feet Storage  
Temperature Range 90°-160° F (Tolerance +/- 2° F)

#### Electrical

120V, 6.3 Amp, 50/60 Hz, .45 kWh (Avg.),  
1535 BTU/hr (Avg.)



SWC243036-G-2B  
2" base (2B) standard for this model

Additional items available for this model (sold separately):

- Mobile Stand (MS2430)
- Mobile Cabinet Stand (WBC2430)
- Roll Out Basket (WB2430)

### SWC182474

#### Overall Dimensions

20.5" D x 24" W x 75" H  
(520.7mm x 609.6mm x 1905mm)  
2 Adjustable Shelves, 1 Fixed Shelf

#### Chamber I.D.

17" D x 20" W x 61" H  
(431.8mm x 508mm x 1549mm)  
12.0 Cubic Feet Storage  
Temperature Range 90°-160° F (Tolerance +/- 3° F)

#### Electrical

120V, 6.5 Amp, 50/60 Hz, .47 kWh (Avg.),  
1604 BTU/hr (Avg.)



SWC182474-G-4B  
4" base (4B) standard for this model

Additional items available for this model (sold separately):

- Mobile Base (MB1824)
- Roll Out Basket (WB1824)

### SWC243074 *Best Seller!*

#### Overall Dimensions

26.5" D x 30" W x 74.5" H  
(673.1mm x 762mm x 1892mm)  
2 Adjustable Shelves, 1 Fixed Shelf

#### Chamber I.D.

23" D x 26" W x 61" H  
(584.2mm x 660.4mm x 1549mm)  
21.1 Cubic Feet Storage  
Temperature Range 90°-160° F (Tolerance +/- 3° F)

#### Electrical

120V, 6.5 Amp, 50/60 Hz, .47 kWh (Avg.),  
1604 BTU/hr (Avg.)



SWC243074-G-4B  
4" base (4B) standard for this model

Additional items available for this model (sold separately):

- Mobile Base (MB2430)
- Roll Out Basket (WB2430)

# Dual Chamber Units

## Warming Cabinets

### DWC183074T

(T = Chambers in Thirds)

#### Overall Dimensions

20.5"D x 30"W x 74.5"H  
(520.7mm x 762mm x 1892mm)  
12.27 Cubic Feet Storage, 3 Adjustable Shelves

#### Upper Chamber I.D.

17"D x 26"W x 15.25"H  
(431.8mm x 660.4mm x 387.3mm)  
3.45 Cubic Feet Storage  
Temperature Range 90°-160° F (Tolerance +/- 1° F)

#### Lower Chamber I.D.

17"D x 26"W x 34.5"H  
(431.8mm x 660.4mm x 876.3mm)  
8.82 Cubic Feet Storage  
Temperature Range 90°-160° F (Tolerance +/- 2° F)

#### Electrical

120V, 8.3 Amp, 50/60 Hz, .82 kWh (Avg.),  
2047 BTU/hr (Avg.)



DWC183074T-G-4B  
4" base (4B) standard for this model

Additional items available for this model (sold separately):

- Mobile Base (MB1830)
- Roll Out Basket (WB1830)

### DWC243074T *Best Seller!*

(T = Chambers in Thirds)

#### Overall Dimensions

26.5"D x 30"W x 74.5"H  
(673.1mm x 762mm x 1892mm)  
17.21 Cubic Feet Storage, 3 Adjustable Shelves

#### Upper Chamber I.D.

23"D x 26"W x 15.25"H  
(584.2mm x 660.4mm x 387.3mm)  
5.28 Cubic Feet Storage  
Temperature Range 90°-160° F (Tolerance +/- 1° F)

#### Lower Chamber I.D.

23"D x 26"W x 34.5"H  
(584.2mm x 660.4mm x 876.3mm)  
11.93 Cubic Feet Storage  
Temperature Range 90°-160° F (Tolerance +/- 3° F)

#### Electrical

120V, 12.5 Amp, 50/60 Hz, .90 kWh (Avg.),  
3071 BTU/hr (Avg.)



DWC243074T-G-4B  
4" base (4B) standard for this model

Additional items available for this model (sold separately):

- Mobile Base (MB2430)
- Roll Out Basket (WB2430)

### DWC243074E

(E = Equal Chambers)

#### Overall Dimensions

26.5"D x 30"W x 74.5"H  
(673.1mm x 762mm x 1892mm)  
17.64 Cubic Feet Storage, 4 Adjustable Shelves

#### Upper Chamber I.D.

23"D x 26"W x 25"H  
(584.2mm x 660.4mm x 635mm)  
8.65 Cubic Feet Storage  
Temperature Range 90°-160° F (Tolerance +/- 2° F)

#### Lower Chamber I.D.

23"D x 26"W x 26"H  
(584.2mm x 660.4mm x 660.4mm)  
8.99 Cubic Feet Storage  
Temperature Range 90°-160° F (Tolerance +/- 2° F)

#### Electrical

120V, 12.5 Amp, 50/60 Hz, .90 kWh (Avg.),  
3071 BTU/hr (Avg.)



DWC243074E-G-4B  
4" base (4B) standard for this model

Additional items available for this model (sold separately):

- Mobile Base (MB2430)
- Roll Out Basket (WB2430)

*Optional*  
**Mobile Cabinet Stand**

Available for 36"H & smaller  
single chamber models

## Contact Us

Contact our customer service department toll free at **877-828-9975**, or by email at [sales@macmedical.com](mailto:sales@macmedical.com).



# Triple Chamber Units

## Warming Cabinets

### TWC183078

#### Overall Dimensions

20.5" D x 30" W x 78.75" H  
(520.7mm x 762mm x 2000mm)  
10.36 Cubic Feet Storage, 1 Adjustable Shelf

#### Upper Chamber I.D.

17" D x 26" W x 11.25" H  
(431.8mm x 660.4mm x 285.7mm)  
2.88 Cubic Feet Storage  
Temperature Range 90°-110° F (Tolerance +/- 1° F)

#### Intermediate Chamber I.D.

17" D x 26" W x 10.75" H  
(431.8mm x 660.4mm x 273mm)  
2.75 Cubic Feet Storage  
Temperature Range 90°-150° F (Tolerance +/- 1° F)

#### Lower Chamber I.D.

17" D x 26" W x 18.5" H  
(431.8mm x 660.4mm x 469.9mm)  
4.73 Cubic Feet Storage  
Temperature Range 90°-160° F (Tolerance +/- 1° F)

#### Electrical

120V, 12.9 Amp, 50/60 Hz, .90 kWh (Avg.),  
3071 BTU/hr (Avg.)



TWC183078-G-4B  
4" base (4B) standard for this model

Additional items available for this model (sold separately):

- Mobile Base (MB1830)
- Roll Out Basket (WB1830)

### TWC243078

#### Overall Dimensions

26.5" D x 30" W x 78.75" H  
(673.1mm x 762mm x 2000mm)  
14.01 Cubic Feet Storage, 1 Adjustable Shelf

#### Upper Chamber I.D.

23" D x 26" W x 11.25" H  
(584.2mm x 660.4mm x 285.7mm)  
3.89 Cubic Feet Storage  
Temperature Range 90°-110° F (Tolerance +/- 1° F)

#### Intermediate Chamber I.D.

23" D x 26" W x 10.75" H  
(584.2mm x 660.4mm x 273mm)  
3.72 Cubic Feet Storage  
Temperature Range 90°-150° F (Tolerance +/- 1° F)

#### Lower Chamber I.D.

23" D x 26" W x 18.5" H  
(584.2mm x 660.4mm x 469.9mm)  
6.4 Cubic Feet Storage  
Temperature Range 90°-160° F (Tolerance +/- 1° F)

#### Electrical

120V, 12.9 Amp, 50/60 Hz, .90 kWh (Avg.),  
3071 BTU/hr (Avg.)



TWC243078-G-4B  
4" base (4B) standard for this model

Additional items available for this model (sold separately):

- Mobile Base (MB2430)
- Roll Out Basket (WB2430)

## Warming Cabinet Accessories



Mobile Stands	
Model #	Overall Dimensions
MS1518	18" D x 21" W x 34" H
MS1824	20" D x 28" W x 34" H
MS1830	20" D x 34" W x 34" H
MS2424	26" D x 26" W x 34" H
MS2430	26" D x 34" W x 34" H


Mobile Cabinet Stands	
Model #	Overall Dimensions
WBC1518	18" D x 21" W x 34" H
WBC1824	20" D x 28" W x 34" H
WBC1830	20" D x 34" W x 34" H
WBC2424	26" D x 26" W x 34" H
WBC2430	26" D x 34" W x 34" H

Mobile Bases	
Model #	Overall Dimensions
MB1824	22" D x 27" W x 6.75" H
MB1830	22" D x 33" W x 6.75" H
MB2424	28" D x 27" W x 6.75" H
MB2430	28" D x 33" W x 6.75" H

Roll Out Baskets	
Model #	Overall Dimensions
WB1518	11.88" D x 10.25" W x 5" H
WB1824	16" D x 16.92" W x 5" H
WB1830	15.88" D x 22.98" W x 5" H
WB2424	21.75" D x 16.98" W x 5" H
WB2430	21.75" D x 23.10" W x 5" H

Note: All warming cabinet mobile stands, mobile cabinet stands, mobile bases, and roll out baskets are sold separately.

### Options

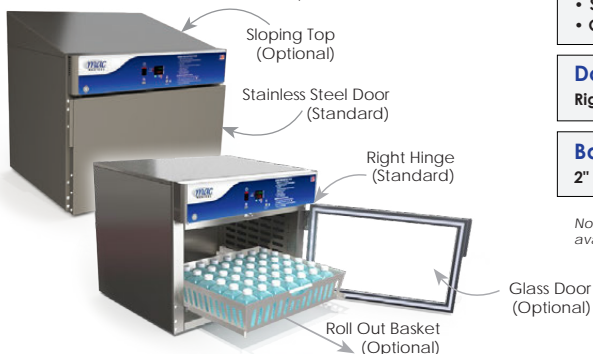

#### Door Types

- Stainless Steel (Standard)
- Glass

#### Door Hinges

- Right (Standard) • Left

#### Base Styles

- 2" • 4" • 6" • No Base

*Note: Additional shelving available for purchase*

#### Options & Accessories


- 220/240V Power Supply
- Celsius
- Two (2) Recessed Unit Options:
  - R1 = Insulation Wrap, No Top Panel, No Side Panels, and No Trim Kit
  - R2 = Top Panel, Side Panels, and Trim Kit
- Seismic Braces
- Mobile Stand
- Mobile Cabinet Stand
- Mobile Base
- Roll Out Basket
- Pass Through Chamber
- Electronic Keypad Lock
- Sloping Top
- Intermediate Chamber Door Lock (only available for Triple Chamber Units)
- Lower Chamber Door Lock (only available for Dual and Triple Chamber Units)

### How to Order Warming Cabinets

Base Model #	Overall Size	Chambers	Door Type	Door Hinge	Base Style	Other Options
SWC151822	17.25"D x 18"W x 22"H	Single	G = Glass Door	LH = Left Hinge	2B = 2" Base	
SWC182424	20.5"D x 24"W x 24.5"H	Single	Blank = Stainless Steel Door (Standard)	Blank = Right Hinge (Standard)	4B = 4" Base	220 = 220/240V Power Supply
SWC243024	26.5"D x 30"W x 24.5"H	Single			C = Celsius	
SWC243036	26.5"D x 30"W x 36.5"H	Single				
SWC182474	20.5"D x 24"W x 75"H	Single				
SWC243074	26.5"D x 30"W x 74.5"H	Single			NB = No Base	R1 = Recessed Unit with Insulation Wrap, No Top Panel, No Side Panels, and No Trim Kit
DWC183074T*	20.5"D x 30"W x 74.5"H	Dual			MB = Mobile Base	
DWC243074T*	26.5"D x 30"W x 74.5"H	Dual				R2 = Recessed Unit with Top Panel, Side Panels, and Trim Kit
DWC243074E**	26.5"D x 30"W x 74.5"H	Dual				SB = Seismic Braces
TWC183078	20.5"D x 30"W x 78.75"H	Triple				WB = Roll Out Basket
TWC243078	26.5"D x 30"W x 78.75"H	Triple				P = Pass Through Chamber
					EL = Electronic Keypad Lock	
					ST = Sloping Top	
					IDL = Intermediate Chamber Door Lock (Only available for Triple Chamber Units)	
					LDL = Lower Chamber Door Lock (Only available for Dual & Triple Chamber Units)	

\* T = Chambers in Thirds

\*\* E = Equal Chambers

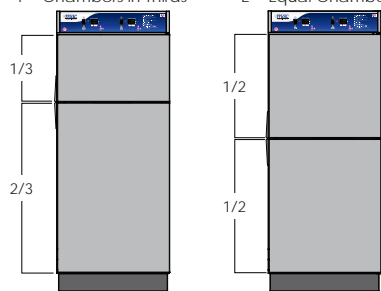

**Example 1:** Model # DWC243074E, Stainless Steel Door, Right Hinge, 4" Base, Celsius, Lower Chamber Door Lock = Order Model # DWC243074E-4B-C-LDL

**Example 2:** Model # SWC182474, Glass Door, Left Hinge, Sloping Top, Mobile Base = Order Model # SWC182474-G-LH-ST-MB

*Note: Model # will change based on features added*





# Transport Stretchers

MAC Medical's Transport Stretchers are designed to ensure safety and efficiency while offering the best blend of features, value, ergonomics, and customer support in the industry. A variety of standard and optional features provide optimal patient comfort & transport.

Best Seller!

# PT1001

## General Transport Stretcher

Economical and versatile, the PT1001 is an ideal solution for any department where cost efficient quality is required. The PT1001 allows for easy and comfortable patient transfer with an increased weight capacity of up to 750 lbs.



**750**  
POUND  
Weight  
Capacity

### STANDARD FEATURES

- 750 lb. weight capacity
- Height adjustment range: 22"–35" (bottom of mattress)
- Ergonomically located oxygen tank holder
- Effortless steering, pivoting, and cornering with 8" casters (1 directional & 3 swiveling)
- 6 brake/steer pedals
- Dual pedestal lift mechanism
- Lift/Trendelenburg pedals on both sides
- Trendelenburg & reverse Trendelenburg +/- 18°
- Dual Fowler assisted lift back rest with 90° adjustability
- Fold forward side rails with built in push handles
- 29" or 31" bed width
- 4" comfort mattress
- 2 hook removable telescopic IV pole
- 4 IV pole sockets available
- Large storage area underneath for patient belongings
- Full perimeter bumper to protect your walls and fixtures
- 4 restraint belt locations


### OPTIONAL FEATURES

- 5th wheel add-on kit
- Fold down push handles
- Knee gatch
- Paper roll dispenser kit \*
- Detachable foot board \*
- Envelope style rail bumper pads \*


\* Additional lead time required (7-10 business days)

Model # will change based on features added


**PT1001**

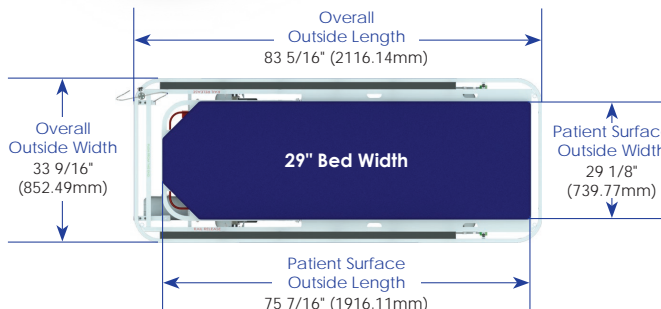

**PT1001-31**

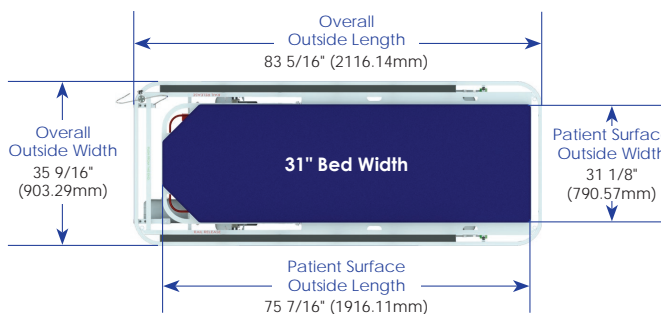

**Height Adjustment Ranges  
(PT1001 & PT1001-31)**

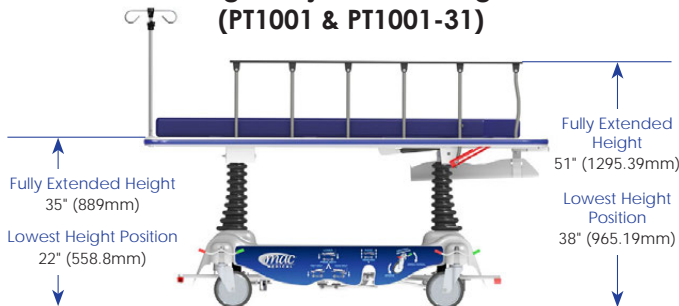

# PTF1000

## 5th Wheel Transport Stretcher

The PTF1000 5th Wheel Transport Stretcher allows for easy and comfortable patient transfer with a weight capacity of up to 750 lbs. The retractable center fifth wheel combines effortless directional movement and cornering with 4 swiveling wheels for ease of use, control, and safety.






### STANDARD FEATURES

- Retractable center fifth wheel provides effortless directional movement and cornering
- 750 lb. weight capacity
- Height adjustment range: 22"-35" (bottom of mattress)
- Ergonomically located oxygen tank holder
- Effortless steering, pivoting, and cornering with 8" casters (4 swiveling)
- 6 brake/steer pedals
- Dual pedestal lift mechanism
- Lift/Trendelenburg pedals on both sides
- Trendelenburg & reverse Trendelenburg +/- 18°
- Dual Fowler assisted lift back rest with 90° adjustability
- Fold forward side rails with built in push handles
- 29" or 31" bed width
- 4" comfort mattress
- 2 hook removable telescopic IV pole
- 4 IV pole sockets available
- Large storage area underneath for patient belongings
- Full perimeter bumper to protect your walls and fixtures
- 4 restraint belt locations

### OPTIONAL FEATURES

- Fold down push handles
- Knee gatch
- Paper roll dispenser kit \*
- Detachable foot board \*
- Envelope style rail bumper pads \*

\* Additional lead time required (7-10 business days)  
Model # will change based on features added

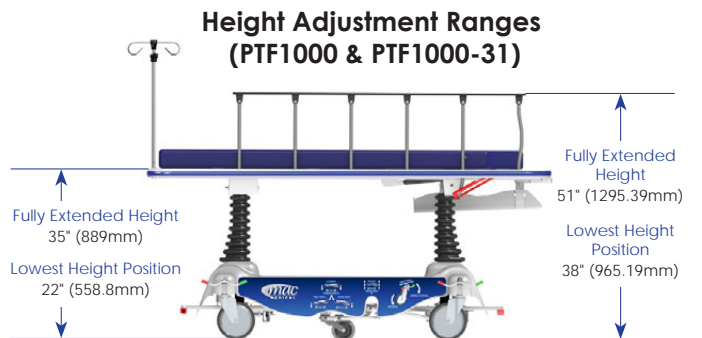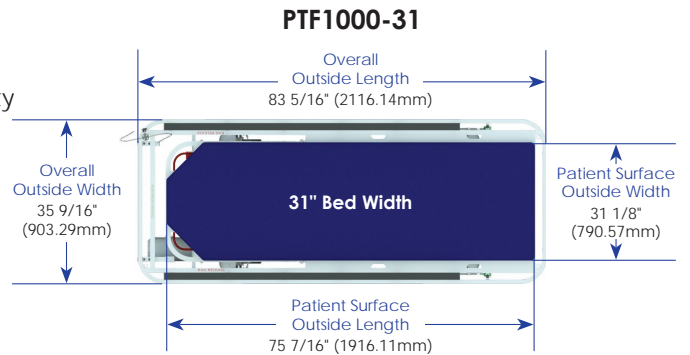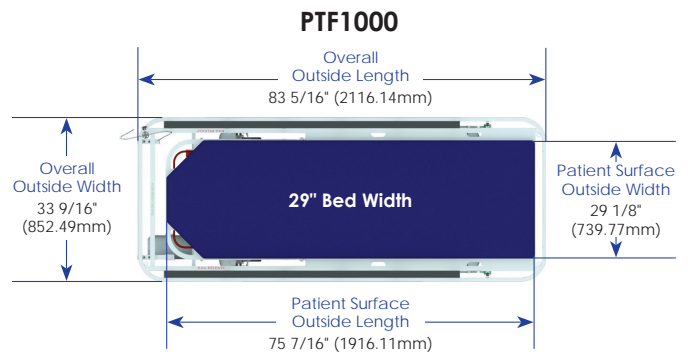

# Which **Stretcher** is right for *your* facility?



Stretcher Comparison Chart		
Standard Features	PT1001	PTF1000
750 lb. weight capacity	X	X
Retractable center fifth wheel		X
Height adjustment range: 22"-35" (bottom of mattress)	X	X
Ergonomically located oxygen tank holder	X	X
8" casters (1 directional & 3 swiveling)	X	
8" casters (4 swiveling)		X
6 brake/steer pedals	X	X
Dual pedestal lift mechanism	X	X
Lift/Trendelenburg pedals on both sides	X	X
Trendelenburg & reverse Trendelenburg +/- 18°	X	X
Dual Fowler assisted lift back rest with 90° adjustability	X	X
Fold forward side rails with built in push handles	X	X
29" or 31" bed width	X	X
4" comfort mattress	X	X
2 hook removable telescopic IV pole	X	X
4 IV pole sockets available	X	X
Large storage area underneath for patient belongings	X	X
Perimeter bumper to protect your walls and fixtures	X	X
4 restraint belt locations	X	X
Optional Features	PT1001	PTF1000
5th wheel add-on kit	X	
Fold down push handles	X	X
Knee gatch	X	X
Paper roll dispenser kit*	X	X
Detachable foot board*	X	X
Envelope style rail bumper pads*	X	X

\* Additional lead time required (7-10 business days)

Model # will change based on features added

## Contact Us

Contact our customer service department toll free at **877-828-9975**, or by email at [sales@macmedical.com](mailto:sales@macmedical.com).



# Surgical Scrub Sinks

Manufactured to the highest quality standards, our surgical scrub sinks are designed to support safe, efficient hand hygiene in healthcare environments. By combining premium materials, reliable components, and expert craftsmanship, each sink is built for durability and reliable functionality in demanding clinical settings.

# ES-Series

## Surgical Scrub Sinks

Our ES-Series Surgical Scrub Sinks were developed to answer space and economical constraints where needed. These sinks have smaller outside dimensions compared to our SS-Series, but still maintain a deep basin to allow for proper scrubbing. The operation is the same as the SS-Series, however features do vary. Choose from single, dual, or triple stations.



### STANDARD FEATURES

- Deep sloping basin
- Adjustable thermostatic controlled pressure regulating mixing valve (one per station, mounted under sink)
- Knee operated water (optional with IR sinks)
- High mount swivel gooseneck spout with adjustable swivel spray aerator
- Flow control valve
- Backflow preventers
- Infrared sinks are equipped with a 24 volt plug in transformer
- Z-brackets (see specifications on rear panel for wall support)
- 304 stainless steel polished and ground to a #4 finish

### OPTIONAL FEATURES

- Infrared operation
- Manual bypass to add knee operated water with IR operation
- Knee operated soap
- Knee operated soap with IR water operation
- Foot operated water
- Foot operated soap
- Digital scrub timer
- Eye wash
- Thermostatic mixing valve (add on to eye wash option)
- Laminar flow
- 220/240V power supply
- Pedestal base

#### ES-Series Standard Models

Base Model #	Overall Dimensions
ES25 (Single Station)	23.5"D x 25"W x 20.5"H
ES47 (Dual Station)	23.5"D x 47"W x 20.5"H
ES76 (Triple Station)	23.5"D x 76"W x 20.5"H

Model # will change based on features added

#### How to Order ES-Series

Infrared Operation	Manual Bypass to add Knee Operated Water with IR Operation	Knee Operated Soap with IR Operation	Foot Operated Water	Foot Operated Soap	Digital Scrub Timer	Eye Wash	Thermostatic Mixing Valve (add on to Eye Wash option)	Laminar Flow	220/240V	Pedestal Base
IR	MB	KSP	FOW	FSP	TM	EYE	TV	LF	220	PB

**Example 1:** Base Model # ES25 with Infrared Operation, Manual Bypass to add Knee Operated Water & Soap, Digital Scrub Timers, and Pedestal Base = Model # **ES25-IR-MB-TM-PB**

**Example 1:** Base Model # ES47 with Laminar Flow = Model # **ES47-LF**


# SS-Series

## Surgical Scrub Sinks


The SS-Series is our premier line of surgical scrub sinks. These sinks combine the best material, components, and craftsmanship available on the market today. Choose from single, dual, or triple stations.



**SS32**  
Single Station

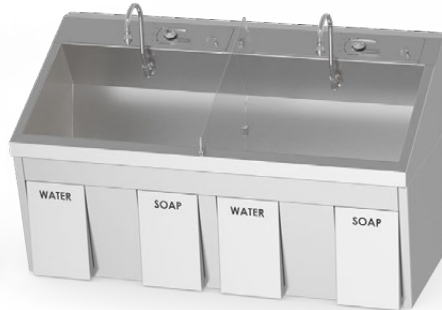

**SS64**  
Dual Station

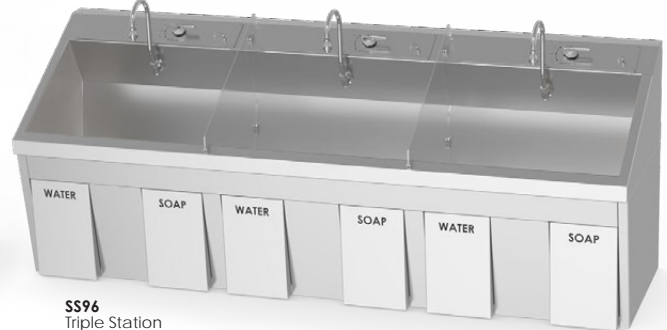

**SS96**  
Triple Station

### STANDARD FEATURES

- Deep sloping basin
- Adjustable thermostatic controlled pressure regulating mixing valve (one per station)
- Knee operated water & soap (optional with IR sinks)
- Check stop strainers for ease of maintenance
- High mount swivel gooseneck spout with adjustable swivel spray aerator
- Flow control valve
- Backflow preventers
- Infrared sinks are equipped with a 24 volt plug in transformer
- Splash screen dividers with dual and triple station sinks
- Main support and Z-brackets (wall support optional)
- 304 stainless steel polished and ground to a #4 finish

### OPTIONAL FEATURES


- Infrared operation
- Infrared soap dispenser
- Manual bypass to add knee operated water & soap with IR operation
- Knee operated soap with IR water operation
- Foot operated water
- Foot operated soap
- Digital scrub timer
- Eye wash
- Thermostatic mixing valve (add on to eye wash option)
- Laminar flow
- 220/240V power supply
- Pedestal base
- In-wall chair carrier
- Wall shelf over sink
- Lighted wall shelf over sink

SS-Series Standard Models	
Base Model #	Overall Dimensions
SS32 (Single Station)	27.5" D x 32" W x 37" H
SS64 (Dual Station)	27.5" D x 64" W x 37" H
SS96 (Triple Station)	27.5" D x 96" W x 37" H

Model # will change based on features added

Recommended Optional Products			
Base Model #	In-Wall Chair Carrier #	Shelf Over Sink #	Lighted Shelf Over Sink #
SS32	S0001	SH0001 (12" D x 32" W)	SH0001-L
SS64	S0002	SH0002 (12" D x 64" W)	SH0002-L
SS96	S0003	SH0003 (12" D x 96" W)	SH0003-L

# SS-Series Features & Options

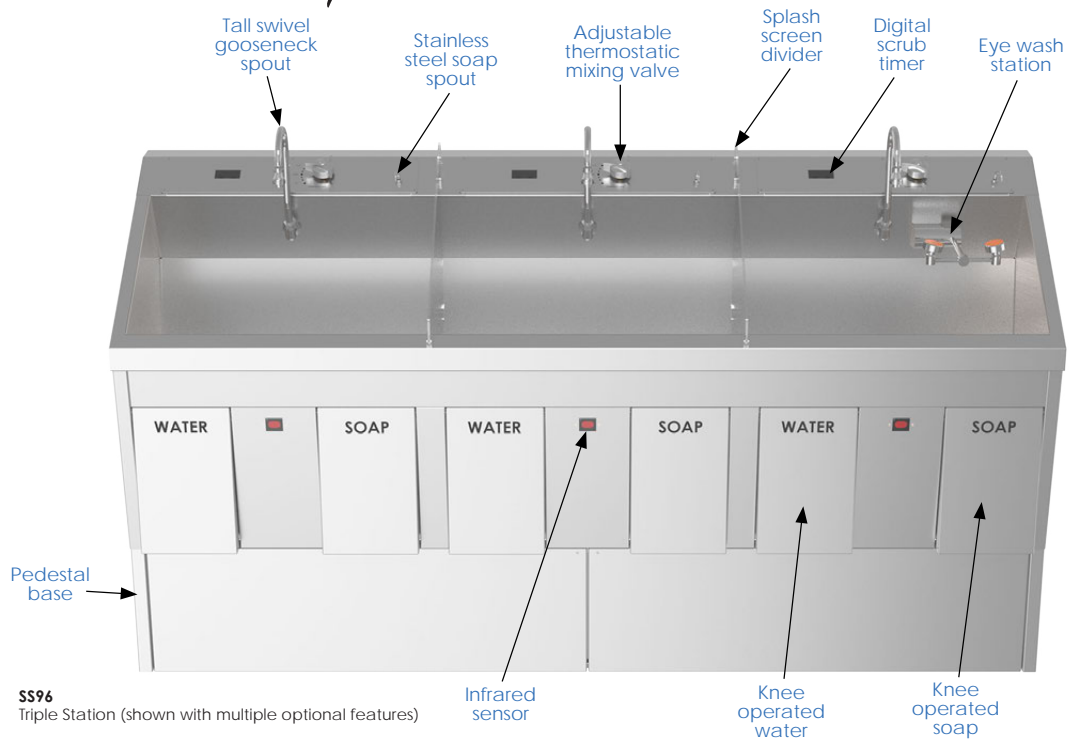

**SS96**  
Triple Station (shown with multiple optional features)

## How to Order SS-Series

Infrared Operation	Infrared Soap Dispenser	Manual Bypass to add Knee Operated Water and Soap with IR Operation	Knee Operated Water with IR Operation	Knee Operated Soap with IR Operation	Foot Operated Water	Foot Operated Soap	Digital Scrub Timer	Eye Wash	Thermostatic Mixing Valve (add on to Eye Wash option)	Laminar Flow	220/240V	Pedestal Base
IR	IRD	MB	KOW	KSP	FOW	FSP	TM	EYE	TV	LF	220	PB

**Example 1:** Base Model # SS32 with Infrared Soap Dispenser, Digital Scrub Timer, and Pedestal Base = Model # **SS32-IRD-TM-PB**

**Example 2:** Base Model # SS96 with Foot Operated Water, Eye Wash, and Laminar Flow = Model # **SS96-FOW-EYE-LF**

## Utility Requirements


Utility	Connection On Unit	Pressure Range	Temperature	Consumption
Drain	1 1/2" x 6" Tail Piece	N/A	N/A	N/A
Cold Water	1/2" NPT	20-50 PSIG	70° F. Max	2 GPM Per Station
Hot Water	1/2" NPT	20-50 PSIG	180° F. Max	2 GPM Per Station
Electrical (IR Units Only)	120V	N/A	N/A	N/A



# Processing Sinks

Designed for superior quality and performance, our processing sinks support efficient cleaning and decontamination in healthcare settings. Crafted with high-quality materials, durable components, and skilled workmanship, each sink provides long-lasting, reliable performance for demanding sterile processing and utility areas.

Ergonomically Friendly

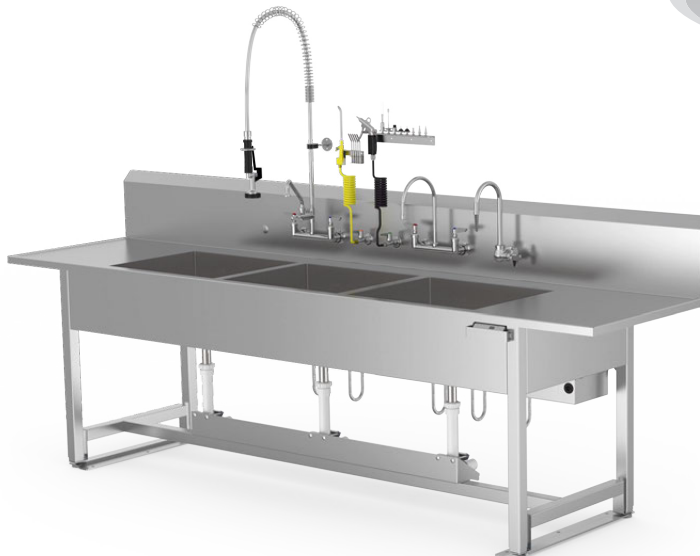

### STANDARD FEATURES

- 12" Angled Backsplash
- 37" Standard Height for Fixed Height Sinks
- Adjustable Feet
- Basin Size: 24"W X 20"L X 10"D
- Etched Gallon Markings
- Flat Working Surface
- Fully Welded Construction with Reinforced Tops
- Marine Edge
- Open Pipe Base
- Seismic Brackets
- Single, Dual, or Triple Basins

**PS3012035-1000**  
 (shown with electronic adjustable height pipe base with universal workflow and optional pre-rinse faucet, mixing valve faucet, DI water faucet, air gun, and water gun)



Note: ETL compliance for adjustable height models only.

### OPTIONAL FEATURES

(Must be specified with sink purchase)

- **Ergonomically friendly** electronic adjustable height hydraulic system with four position programming (electronic adjustable height range: 35"- 42")
- Bolt-on Side Splash
- Copper Drain System (Not Compatible with Deionized Water)
- Hand Lever Drains
- Larger Basin: 30"W X 20"L X 10"D

Electronic Adjustable Height Standard Models	
Base Model #	Overall Dimensions
PS307235-1000	30"D x 72"W x 35"H, Single Basin
PS309635-1000	30"D x 96"W x 35"H, Dual Basin
PS3012035-1000	30"D x 120"W x 35"H, Triple Basin


**Power Requirements (Electronic Adjustable Height Sinks)**  
 120 Volt, 60 Hz, single phase, 15 amp GFCI protected electrical outlet, installed per local building codes, and provides protective earth

Model # will change based on features added

Stationary Height Standard Models	
Base Model #	Overall Dimensions
PS307237-1000	30"D x 72"W x 37"H, Single Basin
PS309637-1000	30"D x 96"W x 37"H, Dual Basin
PS3012037-1000	30"D x 120"W x 37"H, Triple Basin

Model # will change based on features added

## Processing Sink Accessories

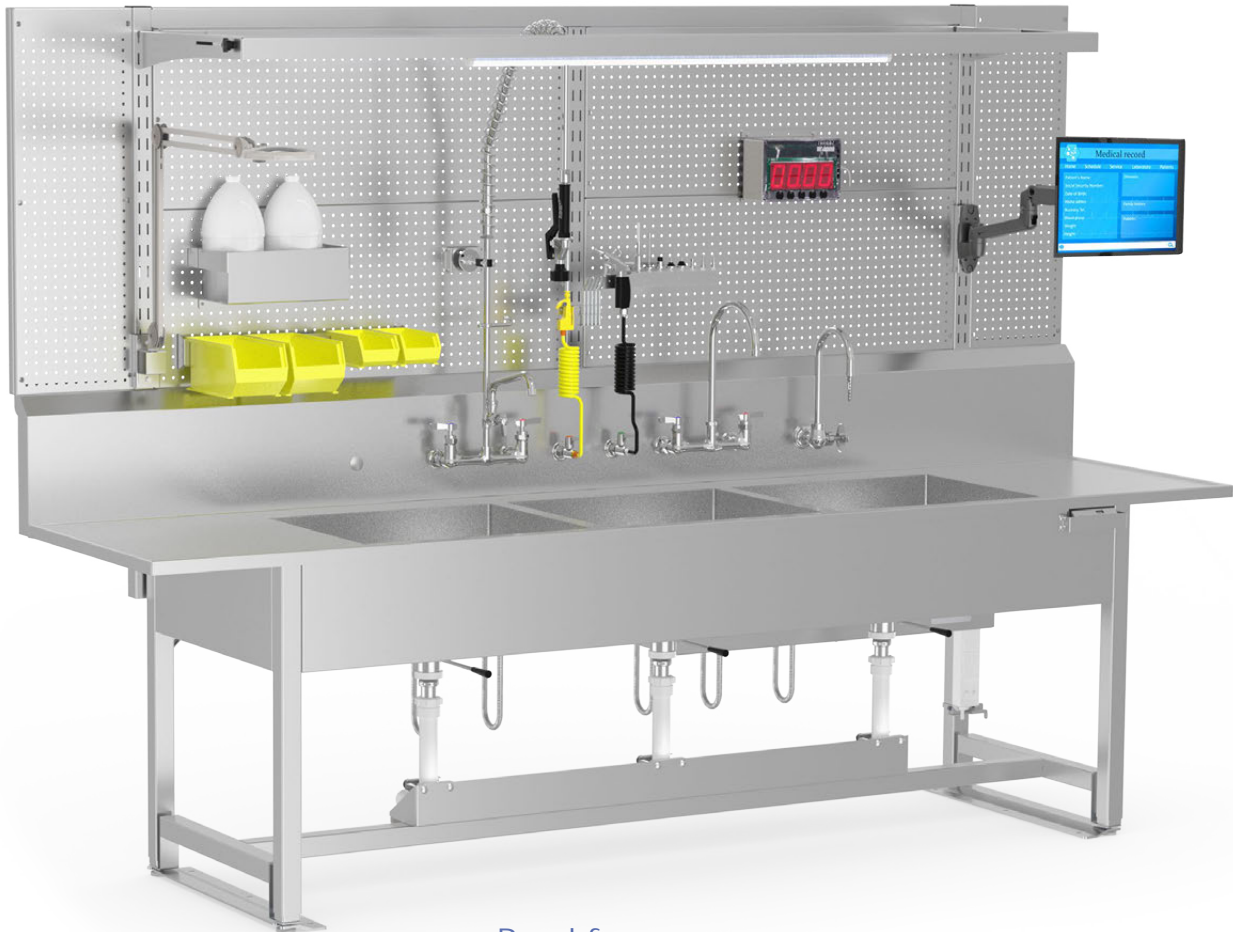

### Backframe w/ optional accessories

**PS3012035-1000**

(shown with electronic adjustable height base with universal workflow and optional pre-rinse faucet, gooseneck faucet, DI water faucet, air gun, water gun, gallon jug holder, LED magnifying inspection lamp, LED soak timer, monitor mount arm, hand lever drains, bin rail with bins, leg mounted power strip, and overhead LED light)

## OPTIONAL SINK ACCESSORIES

- Backsplash Mounted Accessories
  - Air Gun Kit (with 5 Tips)
  - Deionized (DI) Water Faucet
  - Eye Wash
  - Gooseneck Mixing Valve Faucet
  - Pre-Rinse Faucet
  - Water Gun Kit (with 8 Tips)
- Basin Covers
- Leg Mounted Power Strip
- Lower Shelf/Shelves
- Wrist Rest

## PEGBOARD BACKFRAME AND ACCESSORIES

- Pegboard Backframe
  - Bin Rail
    - Plastic Storage Bins – 2 Sizes Available
  - Brush Holder
  - Digital Basin Temperature Indicator
  - Dosing Soap Pump
  - Flushing System
  - Gallon Jug Holder
  - Glove Dispenser
  - Keyboard Arm
  - LED Magnifying Inspection Lamp
  - LED Soak Timer
  - Monitor Mount Arm
  - Overhead LED Light
  - Overhead Light Deflector
  - Pegboard Hooks
  - Pegboard Shelf
  - Spray Bottle Holder
  - Waterproof LED Task Light with Magnifying Lens
  - Work Flow Signs


# Stainless Steel Casework & Cabinets

MAC Medical's full line of stainless steel Wall Cabinets, Base Cabinets, High Cabinets, Operating Room Cabinets, Mobile Cabinets, and Nurses Documentation Stations can be configured in a number of ways to create the arrangement your facility needs.

Whether you need a complete renovation or a new construction, we can work with your elevation plans and facility planners to meet any specific need. As with all our products, all stainless steel cabinets are built with quality and durability in mind.




# From Concept To Completion

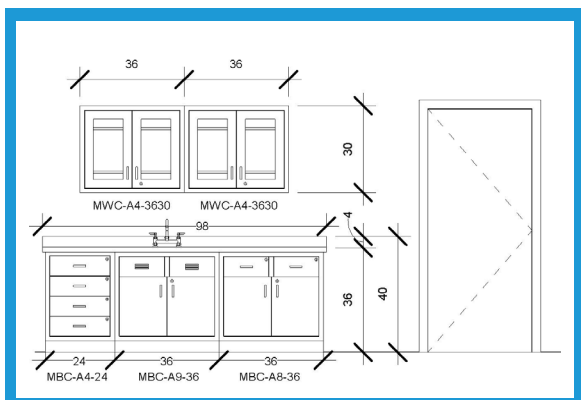

## 1 Plan

MAC Medical's experienced team of engineers can help conceptualize and plan a solution for your casework needs. Whether those needs are simply one cabinet or a complex build-out, we pride ourselves on helping you find a solution that works best for your facility. We will start with a proof of concept, and will revise your drawings until they meet your requirements.

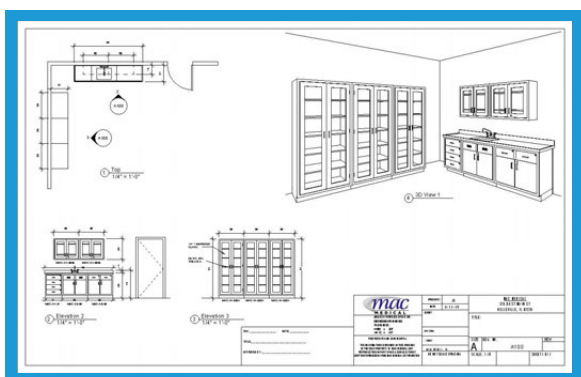

## 2 Build

Once the engineering drawings have been finalized and approved, your casework project starts to come to life in our manufacturing facility. After your project has been built, MAC Medical makes sure your product passes our rigorous quality and standards check list prior to shipment. Quality built and American made continues to be our competitive advantage.



## 3 Install

When your casework arrives at its destination, MAC Medical makes sure installation happens smoothly and timely. We communicate with the contractors/installers to ensure your product is installed according to your specifications. In short, MAC Medical stays involved from concept to completion.

### Contact Us

Contact our customer service department toll free at **877-828-9975**, or by email at [sales@macmedical.com](mailto:sales@macmedical.com).





Wall  
Cabinets

Magnetic door catches  
with a quiet close!



High  
Cabinets



Base  
Cabinets

### STANDARD FEATURES

- Magnetic door catches with a quiet close
- Seamless front face on all cabinets
- Fully welded construction
- Type 300 series stainless steel
- Stainless steel shelves (where applicable, shelves adjust in 1 ½" increments)
- 4" standard pull handle (excludes sliding doors)
- ⅜" tempered glass (where applicable)
- Adjustable feet on Base Cabinets & High Cabinets

*See our casework catalog for model #'s, dimensions,  
and optional features or visit the casework section  
on our website at [www.macmedical.com](http://www.macmedical.com)*



## Operating Room Cabinets

MAC Medical's full line of stainless steel Operating Room Cabinets can be assembled and configured in a number of ways to create a bank of cabinets to fit your facility needs. Standard models are recessed with several optional features to choose from.

### STANDARD FEATURES

- Magnetic door catches with a quiet close
- Seamless front face on all cabinets
- Recessed
- Fully welded construction
- Type 300 series stainless steel
- Stainless steel shelves (excludes Table Accessory Cabinets & Desk Cabinets, shelves adjust in 1 ½" increments)
- ⅜" tempered glass (where applicable)

# Casework **Accessories**

## *Wall Shelves*



### **STANDARD FEATURES**

- Type 300 series stainless steel
- 1" lip on back
- Pre-punched for easy wall mounting (mounting hardware not included)

---

## *Sloping Tops*




### **STANDARD FEATURES**

- Type 300 series stainless steel
- 25° slope angle

## *Corner Guards*



## Counter Tops

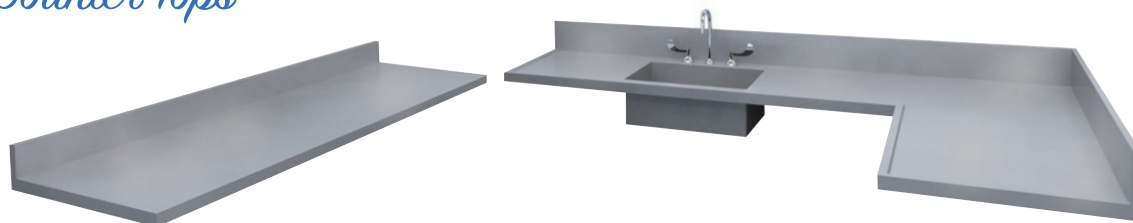

Counter tops are fabricated from type 300 series stainless steel. The underside is reinforced with channels, sink included. Marine edges are provided for counter tops containing sinks. Sinks are fabricated (not drawn) from type 300 series stainless steel with a 5/16" cove corner, laser etched gallon markings, positive pitch towards drain, and are integrally welded to top. A 4" backsplash is standard with optional side splashes. Please consult MAC Medical for custom sizes meeting any and all architectural requirements. Counter tops are fabricated in one piece where design and materials permit.

*All counter tops are sold separately and installed by others along with any plumbing fixtures*

---

## Operating Room Bases - 12" High

### STANDARD FEATURES

- Type 300 series stainless steel

# Nurses Documentation Stations



NDS908430

## STANDARD FEATURES

- Free standing
- Fully welded construction
- Configuration includes wall cabinet, base cabinet, and tower cabinet (tower cabinet can be configured on either side)
- Magnetic door catches
- Stainless steel adjustable shelves
- Standard handles
- Integrated wire management system
- Full CAD renderings
- Engineered room layouts in CAD based on customer provided drawings
- Multiple layouts for individual requirements
- Swivel desk with (2) 5" swivel casters

## OPTIONAL FEATURES

- Stainless steel doors
- Glass hinged doors
- Vented doors for electronics
- No doors
- Pull out writing surface
- Pull out bottom shelf for easy PC install
- Heavy duty pull out drawers
- Key lock
- Tee handles
- Tee handles with key lock
- Electronic keypad door locks
- Pegboard back
- Sloping top (height varies on depth)
- Variable height bases

Standard Model		
Base Model #	Overall Dimensions	Unit Configuration
NDS908430	30"D x 90" W x 84" H	Wall Cabinet (13"D x 48"W x 30"H) Base Cabinet (30"D x 24"W x 32.5"H) Tower Cabinet (30"D x 42"W x 84"H) Counter Top (30"D x 48"W x 36"H) Swivel Desk (22"D x 42"W x 34"H)

*Additional shelving available for purchase, shelving adjusts in 1 1/2" increments*

*Model # will change based on features added*

# Mobile Cabinets

---

## STANDARD FEATURES

- Fully welded construction
- Glass hinged doors
- Stainless steel adjustable shelves
- Standard 4" pull handle
- 4" casters (2 swivel & 2 locking)
- Pull out drawers
- Flat top

## OPTIONAL FEATURES

- Sloping top
- No doors
- Push handles
- Non-marring bumper
- Cam lock
- Tee handles
- Tee handles with key lock
- Electronic keypad door lock
- Pegboard back



Shown with optional sloping top

## Contact Us

Contact our customer service department toll free at **877-828-9975**, or by email at [sales@macmedical.com](mailto:sales@macmedical.com).



# Stainless Steel Carts

MAC Medical manufactures a full line of stainless steel carts that feature fully welded construction for strength and durability. Choose from standard designs or have our engineering department custom design any specific case cart to satisfy your facility needs.

# Closed Case Carts

## STANDARD FEATURES

- Fully welded construction with inner skeletal system for strength and our easy service bolt-on features
- Soundproof tops and doors for quiet operation
- 6" maintenance-free, high temp, high pressure stainless steel casters with sealed ball bearings (3 swivel & 1 fixed directional)
- Bolt-on exterior push handle for better control
- Bolt-on full wrap-around non-marring bumper assembly
- Bolt-on door hinges
- Bolt-on positive lock door latches
- Bolt-on 8.5" x 11" card holder
- Bolt-on tamper evidence security hasps



## OPTIONAL FEATURES

- Glass door(s)
- Etching for serial tag ID/tracking
- Etching for end user logo
- Adjustable solid shelf (non roll out)
- Adjustable solid shelf (roll out)
- Adjustable perforated shelf (non roll out)
- Adjustable perforated shelf (roll out)
- Adjustable wire shelf (non roll out)
- Adjustable wire shelf (roll out)
- Side extensions
- Flag system (one each of red, yellow and green)
- Clean/dirty indicator (slide or wheel)
- Total locking casters
- Tow hitch

### Optional Glass Door Case Cart



### Optional Clean/Dirty Indicator

Available for any new or pre-existing Closed Case Cart  
(Part # CC0230-01-BIO)



## Standard Models

Base Model #	Case Dimensions	Overall Dimensions	Doors
MMCCB3	25" D x 40.25" W x 32" H	26.875" D x 45.75" W x 39.625" H	2
MMCCB4	25" D x 40.25" W x 47" H	26.875" D x 45.75" W x 54.625" H	2
MMCCB6	25" D x 23.625" W x 32" H	26.875" D x 29.125" W x 39.625" H	1
MMCCB9	25" D x 23.625" W x 47" H	26.875" D x 29.125" W x 54.625" H	1
MMCCB12	25" D x 36" W x 32" H	26.875" D x 41.5" W x 39.625" H	2
MMCCB14	25.875" D x 46.5" W x 32" H	27.75" D x 52" W x 39.625" H	2
MMCCB15	25.875" D x 46.5" W x 47" H	27.75" D x 52" W x 54.625" H	2

Model # will change based on features added

# Sealed Case Carts



## STANDARD FEATURES

- Fully welded construction
- Sound-proofed tops and doors for quiet operation
- Exterior push handle for better control
- Door gaskets for increased dust control
- Locking rods to keep doors secure
- Full wrap around non-marring heavy-duty bumper
- Magnetic door catches to keep doors open while in use
- 8.5" x 11" card holder
- Tamper evidence security tabs
- 6" maintenance-free, high-temperature, high-pressure stainless steel casters with sealed ball bearings (2 rigid & 2 swivel with brake)

## OPTIONAL FEATURES

- Etching for serial tag ID/tracking
- Etching for end user logo
- Adjustable solid shelf (non roll out)
- Adjustable solid shelf (roll out)
- Adjustable perforated shelf (non roll out)
- Adjustable perforated shelf (roll out)
- Adjustable wire shelf (non roll out)
- Adjustable wire shelf (roll out)
- Side extensions
- Flag system (one each of red, yellow and green)
- Clean/dirty indicator (slide or wheel)
- Total locking casters
- Tow hitch

Standard Models			
Base Model #	Case Dimensions	Overall Dimensions	Doors
MMSCCB3	25.5" D x 41.25" W x 32" H	27.375" D x 46.75" W x 39.5" H	2
MMSCCB4	25.5" D x 41.25" W x 47" H	27.375" D x 46.75" W x 54.5" H	2
MMSCCB6	25.5" D x 24.625" W x 32" H	27.375" D x 30.125" W x 39.5" H	1
MMSCCB9	25.5" D x 24.625" W x 47" H	27.375" D x 30.125" W x 54.5" H	1
MMSCCB12	25.5" D x 37" W x 32" H	27.375" D x 42.5" W x 39.5" H	2
MMSCCB14	26.375" D x 47.5" W x 32" H	28.25" D x 53" W x 39.5" H	2

*Model # will change based on features added*

**Ergonomically Friendly**

# Vertical Handle Case Carts

## STANDARD FEATURES

- **Ergonomically friendly** vertical handles for precise control
- Fully welded construction
- 8.5" x 11" card holder
- Sound-proofed tops and doors for quiet operation
- Adjustable perforated roll out stainless steel shelves
- Full wrap around non-marring bumper
- Tamper evidence security tabs
- Positive lock door latches
- 6" maintenance-free, high-temperature, high-pressure stainless steel casters with sealed ball bearings (3 swivel & 1 directional lock)

## OPTIONAL FEATURES

- Glass door
- Etching for serial tag ID/tracking
- Etching for end user logo
- Adjustable solid shelf (roll out)
- Adjustable solid shelf (non roll out)
- Adjustable perforated shelf (non roll out)
- Adjustable perforated shelf (roll out)
- Adjustable wire shelf (roll out)
- Adjustable wire shelf (non roll out)
- Side extensions
- Flag system (one each of red, yellow, and green)
- Improved clean/dirty indicator (slide or wheel)
- Tow hitch
- Total locking casters
- Silicon door gaskets for complete sterility control
- Cool touch handles



*Optional*  
**Cool Touch Handles**

Available for any of our Vertical Handle Case Carts (or any 1" diameter vertical handle)

Retrofit kit includes two Cool Touch Handles and corresponding hardware (Part # WC011-100-KT)

*Optional*  
**Clean/Dirty Indicator**

Available for any new or pre-existing Vertical Handle Case Cart (Part # CC0230-01-BIO)

Standard Models		
Base Model #	Case Dimensions	Overall Dimensions
MMCV95	26.5"D x 39.5"W x 40.75"H	27.75"D x 48.25"W x 53"H
MMCV107	24.75"D x 40"W x 38.75"H	26"D x 48.75"W x 51"H

*Model # will change based on features added*

## Open Case Carts

**MMOCC16**  
(shown with 3 optional adjustable perforated shelves)



**MMOCC19**  
(shown with an optional adjustable solid shelf)



### STANDARD FEATURES


- Open design
- Sturdy built cart – perfect for storage, supplies, linens, or distribution
- Full wrap around non-marring bumper
- Stainless steel guard rail
- 6" composite casters\* (2 swivel & 2 locking)

\* Not cart-washable

### OPTIONAL FEATURES

- Adjustable solid shelf (non roll out)
- Adjustable perforated shelf (non roll out)
- Adjustable wire shelf (non roll out)
- Fixed shelving
- Wire grid on back side
- 6" maintenance-free, high temperature, high-pressure stainless steel casters with sealed ball bearings (3 swivel & 1 directional lock)
- Total locking casters
- 6" plastic bin
- 10" plastic bin

**MMOCC21**  
(shown with 2  
optional adjustable  
perforated shelves)

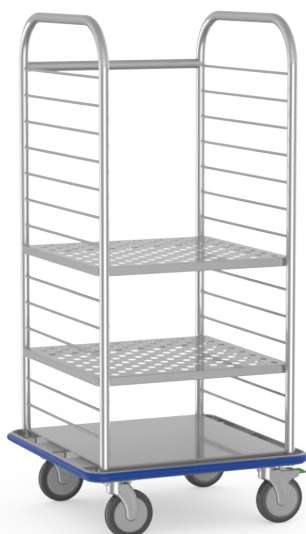

Standard Models	
Base Model #	Overall Dimensions
MMOCC13	29.375"D x 40.375"W x 37.5"H
MMOCC16	26.125"D x 60.375"W x 67"H
MMOCC18	27"D x 40.75"W x 37"H
MMOCC19	30.875"D x 40.75"W x 39"H
MMOCC21	29.75"D x 29.75"W x 62"H
MMOCC24	24"D x 30"W x 54"H
MMOCC25	24"D x 24"W x 54"H

*Shelves ordered separately*

*Model # will change based on features added*

*Some values are rounded up to the next 1/8" for simplicity*

# Utility Carts



## UC *Light Duty*

### STANDARD FEATURES

- Fully welded stainless steel construction
- 3-sided 1" lip (on top & middle shelf) to keep items from rolling off
- Front side turned down (on top & middle shelf) for easier access
- Open design for easy cleaning & sanitizing
- 4" swivel casters (non-locking)
- Push handle
- 200 lb. weight capacity

### OPTIONAL FEATURES

- Locking casters
- Additional shelves
- Customized shelf size & spacing



## AL-200 *Medium Duty*

### STANDARD FEATURES

- Fully welded stainless steel construction
- 3-sided 1" lip (on top & middle shelf) to keep items from rolling off
- Front side turned down (on top & middle shelf) for easier access
- Open design for easy cleaning & sanitizing
- 6" swivel casters (non-locking)
- Push handle
- 300 lb. weight capacity

### OPTIONAL FEATURES

- Locking casters
- Additional shelves
- Customized shelf size & spacing

Standard Models	
Base Model #	Overall Dimensions
UC162432	16"D x 24"W x 32"H
UC162834	16"D x 28"W x 34"H
UC224034	22"D x 40"W x 34"H
AL-200 (A-Frame Design)	19"D x 38.5"W x 35.5"H
AL-200-2442 (A-Frame Design)	24"D x 42"W x 35.5"H

Model # will change based on features added

# Sterile Wrap Carts

## STANDARD FEATURES

- Stainless steel welded construction
- Solid bottom shelf for additional storage
- Four (4) adjustable bars  
(adjustable in 1" increments for comfortable work height)
- 3" casters (2 swivel & 2 locking)



MMSC202863

Standard Model	
Base Model #	Overall Dimensions
MMSC203863	20"D x 38"W x 63"H

*Model # will change based on features added*

### Contact Us

Contact our customer service department toll free at **877-828-9975**, or by email at **sales@macmedical.com**.





◀ Table of Contents

# Stainless Steel Tables

MAC Medical manufactures a full line of stainless steel tables including Instrument Tables, Utility Tables, Nested Tables, Work Tables, Space Saver Tables, Adjustable Over Operating Tables, Central Supply Tables, and Prep & Pack Workstations.

Ergonomically Friendly

# Prep & Pack Workstations

Designed with ergonomics in mind, the Prep & Pack Workstations feature a modular design, allowing for multiple configurations and adjustability for all body types. Choose from the many optional features including adjustable height (manual or electronic), overhead light, articulating magnification light, CPU holder, electrical strip, and more.

## STANDARD FEATURES

- 100% stainless steel construction
- Modular design allows for multiple configurations
- Integrated footrest
- 5" casters (2 swivel & 2 locking)

## OPTIONAL FEATURES

- **Ergonomically friendly** adjustable heights (manual or electronic)
  - Adjustable height range: 34"–46"
- Bin rail
- Overhead light
- Articulating magnification light
- Electrical strip
- CPU holder
- Pull out keyboard tray
- Magnetic dry erase board
- Adjustable flat shelf for riser frame
- Adjustable 45° angled shelf for riser frame
- Lower shelf
- Spool holder
- Drawer boxes (single or double)
- Pegboard panel



**TWS-3672-EAH-B**  
(shown with multiple optional features)

## Standard Models

Base Model #	Overall Dimensions (height at lowest setting)
TWS-3060-B	32"D x 60"W x 84"H
TWS-3072-B	32"D x 72"W x 84"H
TWS-3660-B	38"D x 60"W x 84"H
TWS-3672-B	38"D x 72"W x 84"H

*Model # will change based on features added*

## Utility Tables

### STANDARD FEATURES

- Fully welded construction
- 1 under shelf
- Drawer(s)
- Top guard rail (on 3 sides)
- 3" casters (2 swivel & 2 locking)

### OPTIONAL FEATURES

- Multiple under shelves for one drawer tables
- No guard rail



Standard Models		
Base Model #	Overall Dimensions	Drawers
T0086 <i>Best Seller!</i>	16"D x 20"W x 34"H	1
T0077	16"D x 20"W x 34"H	2

*Model # will change based on features added*

## Nested Tables



### STANDARD FEATURES

- Fully welded construction
- Stainless steel legs (1 1/4")
- Stainless steel tops (type 300 series)
- Stainless steel U brace (1 1/4")
- 2" swivel casters

Standard Models	
Model #	Overall Dimensions
T0089	Set of 6 (T0090 through T0095)
T0090	24"D x 48"W x 42"H
T0091	22"D x 44"W x 40"H
T0092	20"D x 40"W x 38"H
T0093	18"D x 36"W x 36"H
T0094	16"D x 32"W x 34"H
T0095	14"D x 28"W x 32"H

*Model # will change based on features added*

Ergonomically Friendly

## Work Tables



### STANDARD FEATURES

- Fully welded construction
- Reinforced stainless steel top (type 300 series)
- Stainless steel under shelf (type 300 series)
- 2" square legs
- Added supports
- 5" casters (2 swivel & 2 locking)

### OPTIONAL FEATURES

- **Ergonomically friendly** adjustable heights (manual or electronic)
  - Manual adjustable height range: 36"–55"
  - Electronic adjustable height range: 36"–43"
- H or U brace instead of under shelf
- Full base unit underneath
- Drawers under all sizes
- Stainless steel guard rails
- Multiple under shelves
- Adjustable feet

### Standard Models

Base Model #	Overall Dimensions
WT0003	24"D x 36"W x 34"H
WT0004	24"D x 48"W x 34"H
WT0005	24"D x 60"W x 34"H
WT0054	30"D x 72"W x 34"H
WT0055	36"D x 72"W x 34"H
WT0059	30"D x 60"W x 34"H

Model # will change based on features added

Ergonomically Friendly

## Space Saver Tables



### STANDARD FEATURES

- Fold down top shelf with lip (space saving feature!)
- One under shelf
- 2" square legs
- 4" casters (2 swivel & 2 locking)

### OPTIONAL FEATURES

- **Ergonomically friendly** adjustable heights (manual or electronic)
  - Manual adjustable height range: 36"–55"
  - Electronic adjustable height range: 36"–54"

### Standard Models

Base Model #	Overall Dimensions
TSPS-1000-6030	30"D x 60"W x 36"H
TSPS-1000-7230	30"D x 72"W x 36"H

Model # will change based on features added

**Ergonomically Friendly**

## Adjustable Over Operating Tables

### STANDARD FEATURES

- **Ergonomically friendly** manual adjustable height
  - Manual adjustable height range: 36.5"–48.25"
- 3" casters (2 swivel & 2 locking)

### OPTIONAL FEATURES

- 1" lip on back



Standard Models	
Base Model #	Overall Dimensions (height at lowest setting)
T0082-3442	42"D x 34"W x 36.5"H
T0082-4824	24"D x 48"W x 36.5"H
T0082-6024	24"D x 60"W x 36.5"H

*Model # will change based on features added*

## Central Supply Tables

### STANDARD FEATURES

- Heavy duty, stainless steel construction
- Sterile frame
- 2" square legs
- One under shelf
- 5" casters (2 swivel & 2 locking)

### OPTIONAL FEATURES

- 6" plastic bin
- 10" plastic bin
- Drawer



Standard Models	
Base Model #	Overall Dimensions
CST3260-SWF	32"D x 60"W x 34"H
CST2472-SWF	24"D x 72"W x 34"H
CST3684-SWF	36"D x 84"W x 34"H

*Model # will change based on features added*


# Instrument Tables

## STANDARD FEATURES

- Fully welded construction
- Reinforced tops for tables over 25" wide
- Stainless steel legs (1 1/4")
- Stainless steel top (type 300 series)
- Stainless steel under shelf (type 300 series)
- 3" swivel casters

## OPTIONAL FEATURES

- Locking casters
- Guard rails on 3 sides (available on top or bottom shelf)
- Guard rails on all 4 sides (available on top or bottom shelf)
- H or U brace instead of under shelf
- 1" lip on back
- Multiple under shelves

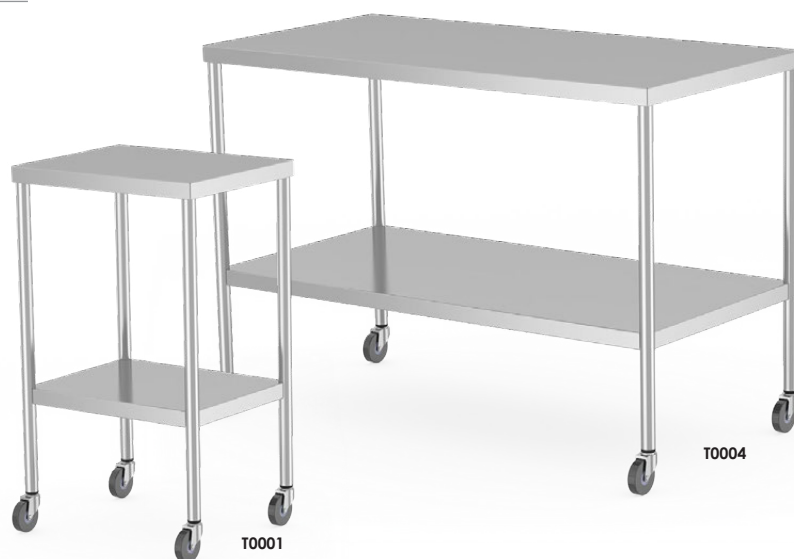

Standard Models	
Base Model #	Overall Dimensions
T0001 <i>Best Seller!</i>	16" D x 20" W x 34" H
T0002	16" D x 30" W x 34" H
T0003	24" D x 36" W x 34" H
T0004 <i>Best Seller!</i>	24" D x 48" W x 34" H
T0005	24" D x 60" W x 34" H
T0006	24" D x 72" W x 34" H

*Model # will change based on features added*

### Contact Us

Contact our customer service department toll free at **877-828-9975**, or by email at [sales@macmedical.com](mailto:sales@macmedical.com).





[« Table of Contents](#)

# Stainless Steel Stands

MAC Medical manufactures a full line of stainless steel stands including Mayo Stands, IV Stands, and Solution Stands.

*Redesigned*

# Mayo Stands

## STANDARD FEATURES

- Full welded stainless steel construction
- Height adjustment range: 36"-58"
- 2" swivel casters



### MYO-2000 *Best Seller!*


- Removable 21 1/4" x 16 1/4" stainless steel tray
- **FOOT** operated friction lock height adjustment/dual post

### MYO-2001

- Same as MYO-2000 except 19 1/8" x 12 5/8" stainless steel tray

### MYO-2002

- Same as MYO-2000 except 25" x 20" stainless steel tray

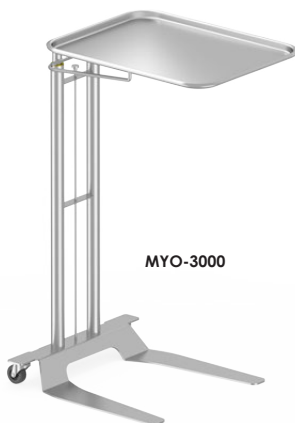

### MYO-3000

- Removable 21 1/4" x 16 1/4" stainless steel tray
- **HAND** operated friction lock height adjustment/dual post

### MYO-3001

- Same as MYO-3000 except 19 1/8" x 12 5/8" stainless steel tray

### MYO-3002

- Same as MYO-3000 except 19 1/8" x 12 5/8" stainless steel tray



# IV Stands

## STANDARD FEATURES

- Stainless steel construction
- Height adjustment range: 68"–108"
- Tip resistant base
- 2" swivel casters



### IVS-3000 Series

- Foot operated
- 5-leg stand
- Welded base

#### IVS-3004

- Four hook design

#### IVS-3006

- Six hook design



### IVS-5000 Series


- Hand operated
- 4-leg stand
- Knock down version (assembly required)

#### IVS-5004

- Four hook design

#### IVS-5006

- Six hook design

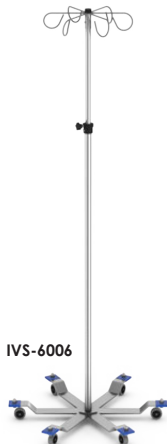

### IVS-6000 Series

- Hand operated
- 6-leg stand
- Knock down version (assembly required)

#### IVS-6004

- Four hook design

#### IVS-6006

- Six hook design



### IVS-7000 Series

- Hand operated
- 5-leg stand
- Knock down version (assembly required)

#### IVS-7004 *Best Seller!*

- Four hook design

#### IVS-7006

- Six hook design

# Solution Stands

## STANDARD FEATURES

- Fully welded stainless steel construction
- Single and double basin designs available
- Available with lower shelf
- Open frame design allows accessibility from any side
- 2" swivel casters



SOL-1000

### SOL-1000

- Overall Dimensions: 14"D x 14"W x 34"H
- Removable 7 quart stainless steel basin



SOL-1001

### SOL-1001 *Best Seller!*

- Same as SOL-1000 except with stainless steel lower shelf
- Removable 7 quart stainless steel basin



SOL-2000

### SOL-2000

- Overall Dimensions: 16"D x 27.5"W x 34"H
- Two removable 7 quart stainless steel basins



SOL-2001

### SOL-2001

- Same as SOL-2000 except with stainless steel lower shelf
- Two removable 7 quart stainless steel basin

## Contact Us

Contact our customer service department toll free at **877-828-9975**, or by email at [sales@macmedical.com](mailto:sales@macmedical.com).



## Other Stainless Steel Products

MAC Medical manufactures several other stainless steel products including Pass Through Windows, Linen Hampers, Gowning Benches, Step Stools, Kick Buckets, Stationary Wall Shelves, Fold Down Wall Shelves, and Peg Boards.

# Pass Through Windows

## STANDARD FEATURES

- 3/16" tempered glass
- Adjustable width for wall thickness (between 4.5"–6.125")
- 4-sided gasket
- Dual swing hinge with adjustable spring tension

## OPTIONAL FEATURES

- Decontamination side shelf (various sizes)
- Clean side shelf (various sizes)



Standard Models	
Base Model #	Outside Dimensions
PTW0001	37.5"W x 34.5"H (no shelf)
PTW0002	37.5"W x 34.5"H (shelf included)
PTW0003	24"W x 34.5"H (no shelf)

*Model # will change based on features added*

# Kick Buckets

## STANDARD FEATURES

- Stainless steel frame
- Perimeter bumper to protect walls
- Tip resistant
- 2" swivel casters



### KB-1000 *Best Seller!*


- Overall dimensions: 14"H x 13 1/2" diameter
- Removable 13 quart stainless steel basin



### KB-2000

- Overall dimensions: 8"H x 16 1/4" diameter
- Removable 9 1/2 quart stainless steel basin

## Gowning Benches

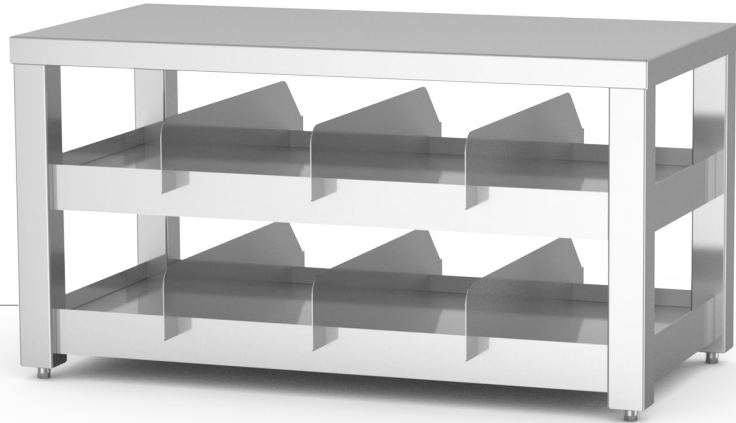

### STANDARD FEATURES

- Stainless steel construction
- H-Brace
- Adjustable feet

### OPTIONAL FEATURES

- Casters
- Dividers
- Middle & bottom shelf

### Standard Models

Base Model #	Overall Dimensions
GB183619-2SH-AF	18" D x 36" W x 19" H
GB184819-2SH-AF	18" D x 48" W x 19" H

*Model # will change based on features added*

## Step Stools



### STANDARD FEATURES

- Heavy duty, stainless steel construction
- Tip resistant
- Corrugated, black rubber top surface on each stool
- Rubber floor tips
- Single step stools (STS121805) are stackable and connect together to create a large platform

### OPTIONAL FEATURES

- Easy roll around transport cart with 4 swiveling casters (2 locking) available for STS121805 – holds up to 5 single step stools

### Optional Transport Cart

Available for STS121805  
(Part # T0800)

Holds up to 5 single step stools



### Standard Models

Model #	Overall Dimensions
STS121805 (single step) <b>Best Seller!</b>	12" D x 18" W x 6.25" H
STS121808-HD (single step with handle)	12" D x 18" W x 6.25" H (41" H at handle)
STS263015 (two steps)	26" D x 30" W x 15" H

# Linen Hampers

## STANDARD FEATURES

- Fully welded, heavy gauge all stainless steel construction
- Heavy gauge platform to support bags and ease of removal
- 2" swivel casters


## OPTIONAL FEATURES

- Foot operated stainless steel lid
- Rotating bumper guards over casters




### LH-1000

- Overall dimensions: 35"H x 57 1/4" perimeter
- Accommodates 18" diameter bags

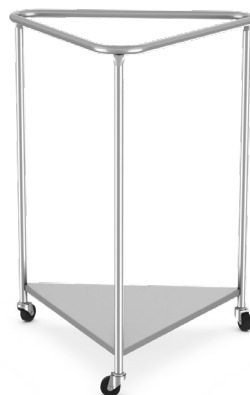

### LH-1001


- Overall dimensions: 35"H x 78 5/8" perimeter
- Accommodates 25" diameter bags

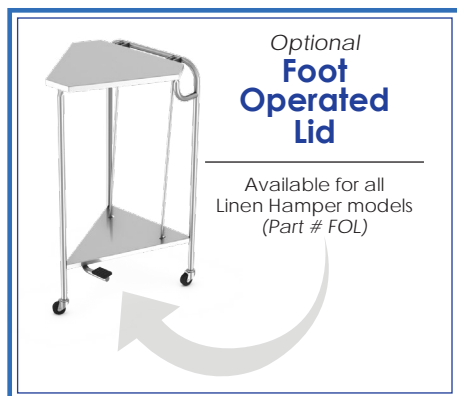

### Optional Foot Operated Lid

Available for all  
Linen Hamper models  
(Part # FOL)

## Stationary Wall Shelves

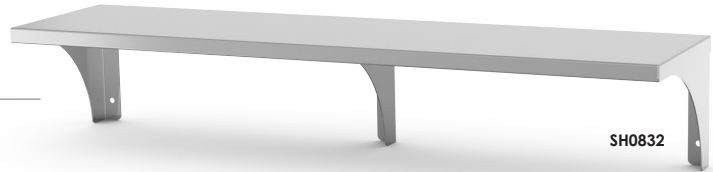

### STANDARD FEATURES


- Stainless steel construction
- Pre-drilled for easy wall mounting (mounting hardware not included)

### OPTIONAL FEATURES

- Bullnose front

Standard Models		
Model #	Overall Dimensions	Brackets
SH0624	6"D x 24"W	2
SH0630	6"D x 30"W	2
SH0824	8"D x 24"W	2
SH0830	8"D x 30"W	2
SH0832	8"D x 32"W	3
SH0836	8"D x 36"W	3
SH1224	12"D x 24"W	2
SH1236	12"D x 36"W	3
SH1248	12"D x 48"W	3
SH1260	12"D x 60"W	4

## Fold Down Wall Shelves

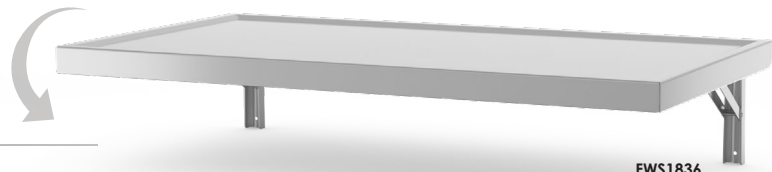

### STANDARD FEATURES

- Stainless steel construction
- Fold down shelf
- Marine edge
- Pre-drilled for easy wall mounting (mounting hardware not included)

### OPTIONAL FEATURES


- Bullnose front

Standard Models		
Model #	Overall Dimensions	Brackets
FWS1224	12"D x 24"W	2
FWS1236	12"D x 36"W	2
FWS1248	12"D x 48"W	2
FWS1824	18"D x 24"W	2
FWS1836	18"D x 36"W	2
FWS1848	18"D x 48"W	2

# Peg Boards

## STANDARD FEATURES


- Stainless steel construction (3/4" overall thickness from wall)
- 9/32" diameter holes on 1" center
- Key hole slots for mounting (mounting hardware not included)
- Can be mounted vertically or horizontally

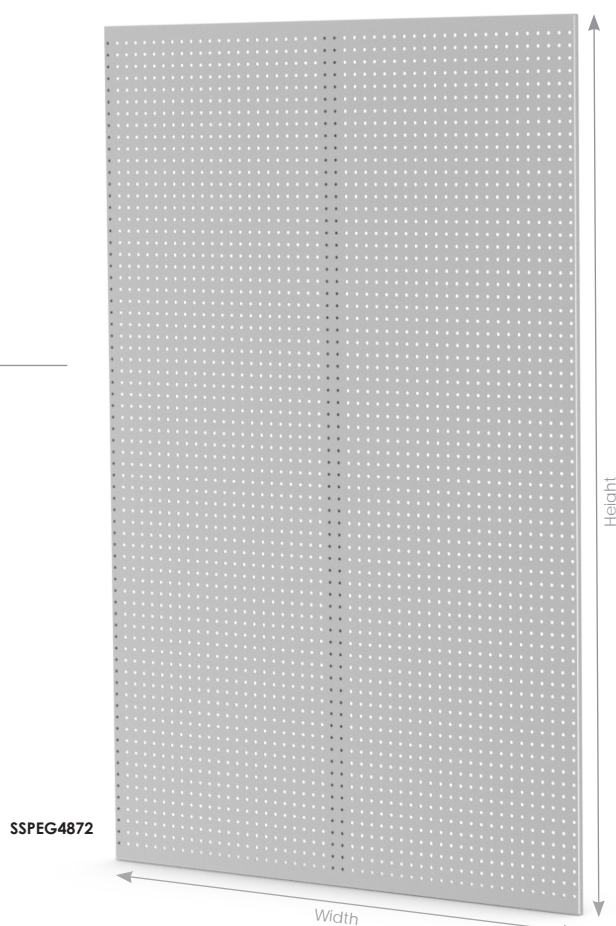

Standard Models	
Model #	Overall Dimensions
SSPEG2436	24"W x 36"H
SSPEG3636	36"W x 36"H
SSPEG3648	36"W x 48"H
SSPEG4248	42"W x 48"H
SSPEG4848	48"W x 48"H
SSPEG4872	48"W x 72"H
SSPEG4896	48"W x 96"H

Standard Hooks	
Model #	Hook Size
IC0128-01	3-3/32" stainless steel single peg board hook (fits 3/16" & 1/4" diameter holes)
IC0128-02	5-3/16" stainless steel single peg board hook (fits 3/16" & 1/4" diameter holes)
IC0128-03	8-3/4" stainless steel single peg board hook (fits 3/16" & 1/4" diameter holes)
IC0901	8" stainless steel single peg board hook (fits 1/4" diameter holes)
IC0128-05	9" stainless steel single peg board hook (fits 1/4" diameter holes)

All hooks are compatible with all peg board models and are sold separately

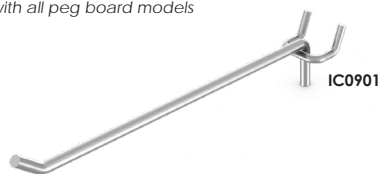



## Laminate Mobile **Medical Carts**

MAC Medical's full line of laminate mobile medical carts can be configured in a number of ways to create the arrangement your facility needs. Several optional features are available also. As with all our products, our laminate mobile medical carts are built with quality and durability in mind.


## STANDARD FEATURES

- 1 1/4" Base with full wrap around non-marring bumper
- Storage panels for organizing trays and baskets diagonally or horizontally
- 4" composite casters (2 swivel & 2 locking)
- Magnetic catches on solid and glass doors
- Sturdy built cart – perfect for storage, supplies, linens, or distribution

## OPTIONAL FEATURES

- Locks (cam locks or digital e-locks)
- Solid surface tops on short laminate mobile carts
- Total locking casters

## MOBILE CART MATERIALS

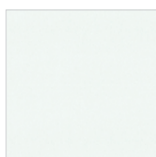

**White TFL Laminate**

Durable thermal-fused white laminate is on 3/4" thick high-density engineered wood



**Gray TFL Laminate**

Durable thermal-fused gray laminate is on 3/4" thick high-density engineered wood



**Aluminum HPL Laminate**

Brushed aluminum high-pressure laminate is a cost effective way to get the appearance of stainless steel at a lower price

## MOBILE CART MODELS

### *Short Mobile Carts*



**MMLCSW245836**

Picture above is a wide gray laminate cart with bow pulls, standard top, a 5th wheel for stability, and optional tambour doors & cam locks



**MMLCSM314036**

Picture above is a medium width gray laminate cart with standard top, wire pulls, and optional glass doors & storage bins



**MMLCSN312436**

Picture above is a narrow white laminate cart with a standard solid door & top, and a wire pull

Short Models (36" Tall)		
Base Model #	Overall Dimensions	Cabinet Dimensions
MMLCSW245836	23 3/4" D x 57 7/8" W x 36" H	18 5/8" D x 53 7/8" W x 29 7/8" H
MMLCSM314036	31 1/8" D x 40 1/2" W x 36" H	26 1/2" D x 36 1/2" W x 29 7/8" H
MMLCSN312436	31 1/8" D x 23 3/8" W x 36" H	26 1/2" D x 18" W x 29 7/8" H

Tall Models (78" Tall)		
Base Model #	Overall Dimensions	Cabinet Dimensions
MMLCTM314078	31 1/8" D x 40 1/2" W x 78" H	26 1/2" D x 36 1/2" W x 71 7/8" H
MMLCTN312478	31 1/8" D x 23 3/8" W x 78" H	26 1/2" D x 18" W x 71 7/8" H

### Tall Mobile Carts



**MMLCTM314078**

Picture above is a medium width gray laminate cart with glass doors, wire pulls, and a digital e-lock



**MMLCTN312478**

Picture above is a narrow aluminum laminate cart with tambour door, bow pull, and a cam lock

### MOBILE CART DOOR STYLES



**Solid**

Solid doors provide the best security. They are strong and more tamper resistant.



**Glass**

Glass doors give you the accessibility of seeing what is in your cart without opening the doors.



**Tambour**

Tambour doors slide open and save space as the doors roll up. You don't have to allow space for the doors to swing open.

## Door Accessories

### Door Pulls



Wire Pull



Flat Pull



Bow Pull

### Locks



Cam Lock



Digital E-Lock

Not available with Tambour doors.

### Magnetic Catches



Glass and solid doors come with magnetic catches.

## MOBILE CART INTERIOR ACCESSORIES

### Trays & Baskets

#### 40x60 Clear Models



**5cm TRAY - MMT054060**  
SHORT DIVIDER - MMDO540  
LONG DIVIDER - MMDO560



**10cm TRAY - MMT104060**  
SHORT DIVIDER - MMD1040  
LONG DIVIDER - MMD1060



**20cm TRAY - MMT204060**  
SHORT DIVIDER - MMD2040  
LONG DIVIDER - MMD2060




**10cm TRAY - MMT104060ND**  
(NOT DIVIDABLE)  
\* LID - L4060



**10cm BASKET - MMB104060**  
SHORT DIVIDER - MMD1040  
LONG DIVIDER - MMD1060



**20cm BASKET - MMB204060**  
SHORT DIVIDER - MMD2040  
LONG DIVIDER - MMD2060

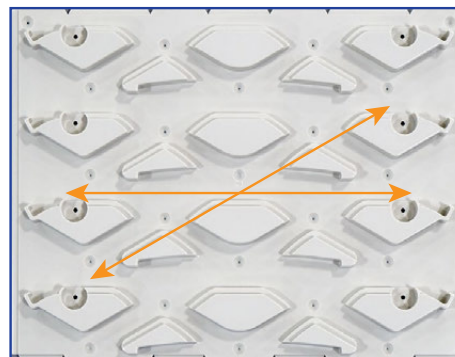

Trays and Baskets can be stored both diagonally and horizontally.



## Company-wide Green Initiative

MAC Medical is committed to making sound decisions to manage our environmental footprint and protect the environment. Our company-wide green initiative is conscious of our impact on the planet and underscores our commitment to sustainability.



### Recycled Material

MAC Medical purchases nearly 1.7 million lbs. of recycled steel each year. Carbon emissions are reduced by 70% when using recycled steel versus newly mined steel. This equates our usage of recycled steel to an approximate 1,884 metric ton reduction in carbon emissions. It is estimated that planting 5 trees reduces 1 metric ton of CO<sub>2</sub>. Thus, our usage of recycled steel versus newly mined steel equates to us planting 9,444 trees last year.

Additionally, our corrugated packaging products have a recycled content of 36%. This is comprised of 30% post-consumer fiber (recycled boxes) and 6% pre-consumer fiber (waste created at the mills that is recycled).



### Facility Lighting

MAC Medical's manufacturing facility has been outfitted with skylights and energy efficient LED lights. Skylights help reduce the amount of electricity needed to illuminate the working space, while LED lights provide a longer operational life span and an overall reduction in wattage used. In fact, 24,298 watts are saved in our entire facility for every hour LED lights are in use. This reduction in wattage reduces carbon emissions and equates to us planting approximately 150 trees per year.



### Manufacturing Practices

MAC Medical practices lean manufacturing to aid in the manufacturing process. Our lean techniques enable us to deliver high quality products at lower costs and ultimately minimizing our impact on the planet. One of those techniques is our investment in fiber laser cutting technology which allows for the maximization of each stainless steel sheet used in production. Even though 100% of our scrap steel is recycled, our fiber laser cutters largely reduce the amount of scrap sent to the smelter.



# Warranty Information

## Limited Lifetime Warranty

MAC Medical warrants to the original purchaser that its Proprietary Products will be free from defects in workmanship or materials under normal use and service for the life of the product, so long as owned by the original purchaser, according to the limitations set forth below. Defective products shall be repaired or replaced at MAC Medical's option at no cost to the original purchaser provided:

1. The customer must obtain a written return authorization supplied by MAC Medical's customer service department. The product must be returned, properly packaged, with a copy of the original sales receipt and copy of the RMA authorization. Please call (877) 828-9975 or (618) 476-3550 to receive a return authorization.
2. The customer must pre-pay freight charges to and from MAC Medical and must ship merchandise properly packaged, in a way the product will not be damaged during transit to the factory. MAC Medical does not warranty any freight damage to or from the factory. It is the customer's responsibility to inspect the product for packaging damage before signing the BOL.

This warranty does not apply to products, which have been subject to abuse, misuse, accident, modification, alteration, tampering, negligence, lack of routine maintenance or misapplication; or products that have been repaired by other than MAC Medical or its authorized representatives.

This warranty does not apply to glass, fabrics, vinyl, seat coverings, cushions, padding or their stitching, gluing or installation. Component parts not manufactured by MAC Medical, such as casters, caster inserts, any components made from rubber or plastic, circuit boards, fan heaters, plumbing parts, electrical switches and other components are also excluded. For these component parts, MAC Medical will pass on the original manufacturer's warranty to MAC Medical equipment original purchasers. Please consult factory for questions regarding the warranties of these component parts. This warranty does not apply to custom fabricated products. Consult factory.

UNLESS EXPRESSLY SET FORTH IN THIS WARRANTY, THERE ARE NO OTHER WARRANTIES, WHETHER EXPRESS OR IMPLIED, OF FITNESS AND / OR MERCHANTABILITY OR ANY OTHER WARRANTY IMPLIED BY CUSTOM, USAGE OR COURSE OF DEALING.

Liability of MAC Medical under this warranty is limited to the repair and / or replacement of any products. MAC Medical specifically excludes and disclaims any responsibility for any incidental or consequential damages claimed to have arisen from any defect in workmanship or materials. No representative has any authority to change or enlarge the above warranty or MAC Medical liability. This supersedes all previous warranty provisions.

# Contact Us

---


## Corporate Office

609 S. Breese St.  
Millstadt, IL 62260

## Manufacturing Plant

820 South Mulberry Street  
Millstadt, IL 62260

Both facilities are located in Southern Illinois, making this an ideal place for centralized distribution to the entire United States.

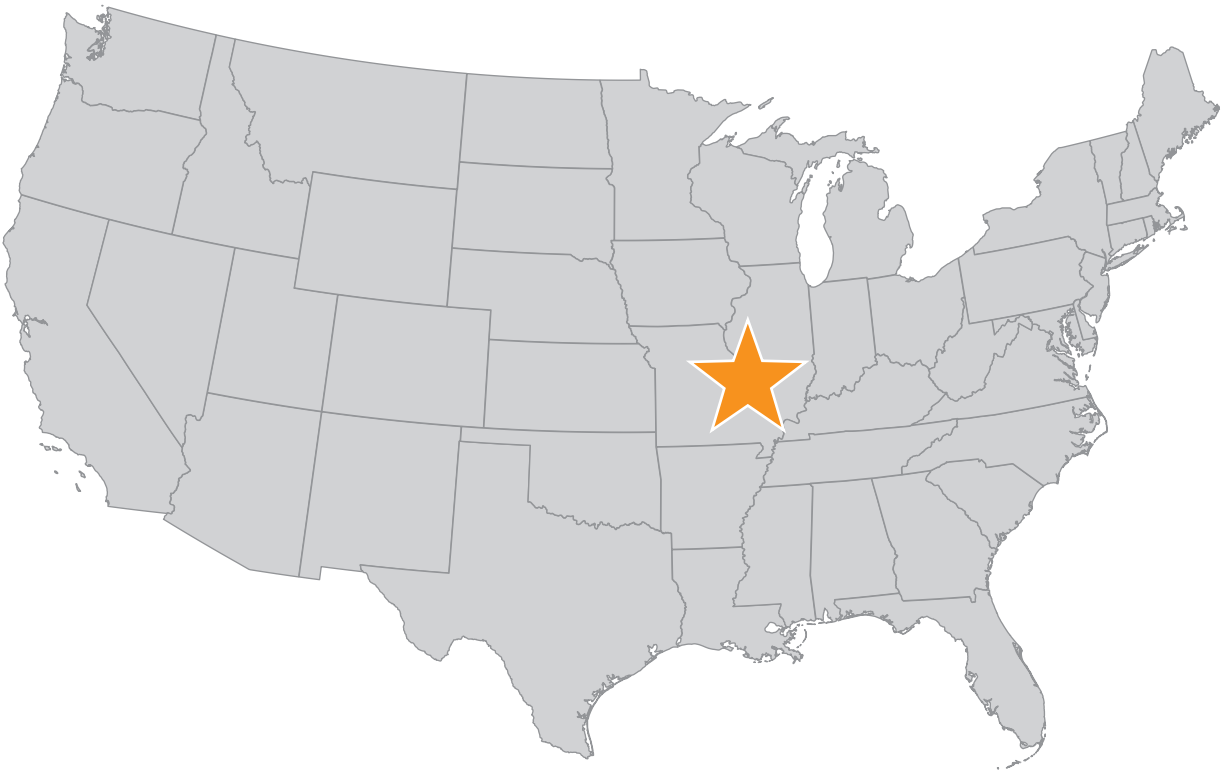

*For questions or to place an order, contact our customer service department:*

**Phone:** 618-476-3550

**Toll Free:** 877-828-9975

**Fax:** 618-476-3337

**Email:** [sales@macmedical.com](mailto:sales@macmedical.com)

**Website:** [www.macmedical.com](http://www.macmedical.com)

## Thank You!

Our company has been in business for 25+ years and we could not have done it without people like you! We are grateful for your partnership and support. From humble beginnings to our current day manufacturing technology, serving the medical community has never been a task for MAC Medical, but more of a pleasure. As our company expands and continues to grow we will strive for excellence in our goal of serving you. Thank you for being a part of our story.



Sincerely,



Dennis Cooper  
CEO/President

INSPIRED  
TO  
BE BETTER

MADE **USA**  
IN THE



MAC Medical, Inc.  
Publication No. CFU-01 Rev H  
March 2026

Information regarding products are subject to change without prior notice.

**Corporate Office**  
609 S. Breese St.  
Millstadt, Illinois 62260

**Manufacturing Plant**  
820 South Mulberry Street  
Millstadt, Illinois 62260

Phone: 618-476-3550 • Toll Free: 877-828-9975 • Fax: 618-476-3337  
Email: [sales@macmedical.com](mailto:sales@macmedical.com) • Website: [www.macmedical.com](http://www.macmedical.com)